Protac MyBaSe®

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GB - Directions for use

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SE - Bruksanvisning

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Protac MyBaSe®

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GB - Directions for use

Protac MyBaSe®

You have selected a quality Protac product. Protac guarantees that you have purchased an original, Danish designed product developed in consultation with experts.

Protac has developed, produced and sold sensory-stimulating aids to the healthcare sector in Denmark, the rest of Europe and Australia since 1994.

Protac products are developed on the basis of theories on sensory integration, and are supported by extensive research, documentation and case studies.

Protac focuses strongly on protection of the environment and safety during its production process. Protac's plastic balls are TÜV-approved and all textiles are OEKO-certified.

GB

For more details, go to www.protac.dk

Labelling

Our products are labelled according to legislation concerning medical equipment.

The label with details of manufacturer, raw materials and washing instructions is sewn inside the cover of the air compartment to avoid irritating the user when lying on the mattress. This location also protects the label from wear and tear.

Key to label symbols: please refer to **Key to symbols**

General precautions - Safety

- ⚠ • **Not for unsupervised use** by users unable to care for themselves.
- ⚠ • **Not for unsupervised use by** users with a poor sense of balance.
- ⚠ • Always carefully monitor a user's reaction.
- ⚠ • **Do not** allow the user to lie with their head down on the mattress.
- ⚠ • **Do not** allow more than 1 user at a time to use the mattress unsupervised.
- ⚠ • The Protac MyBaSe® **is not** designed for outdoor use.
- ⚠ • **Not** to be used in water.
- ⚠ • **Must not** be used near an open flame.
- ⚠ • The Protac incontinence cover should be pulled over the mattress in the event of incontinence problems (urination problems) the same risk **warnings** apply as for the use of plastic bags:
 - mattress with incontinence cover **must not** be used for children without supervision
 - mattress with incontinence cover **must not** be used for people who cannot take care of themselves
 - you **must not** lie with the head downwards.

Product purpose and areas of use

Sensory integration is the idea behind the development of all Protac products. The theory originated in the 1970s in the USA, and was developed by A. Jean Ayres, an occupational therapist and psychologist. The theory is based on how the brain processes sensory impressions, and how the various sensory stimuli can be influenced via applying weight and pressure to the tactile sense (the sense of touch) in the skin and the proprioceptive sense (the sense of movement) in muscles and joints. The pressure and movement of the balls act on these senses in particular.

The Protac MyBaSe® was developed to help sufferers of function impairment.

Children, young and old people can benefit from the mattress:

- Brain damage, including spasticity, Huntington's Chorea, Parkinson's Disease and other ailments
- ADHD and autism
- Psychiatric diseases, such as dementia, schizophrenia, depression, anxiety, anorexia etc.
- Development problems, sensory problems, deaf/blind, etc.
- Stress

The Protac MyBaSe® is used for lying on as an active surface for sufferers of brain damage with spasticity. The pressure of the balls helps suppress spasticity, and calm the body of patients with brain damage. Furthermore, the patient benefits from enhanced body awareness.

The Protac MyBaSe® is used for ADHD and autism sufferers to lie, roll and jump on.

The active side of the mattress stimulates the sense of balance, whilst the pressure of the balls creates calm and better body awareness.

GB

Psychiatric patients often suffer from muscular and psychological unrest. Protac MyBaSe® is used to create calm and comfort for the user when lying on the mattress.

It is used for sufferers of a range of development problems, sensory problems, deaf/blindness and to create enhanced sensory stimuli in the body.

The Protac MyBaSe® can be used in institutions, sheltered housing, day centres, hospitals, patient lounges and private homes. The mattress can be used in the gymnastics hall of schools, for instance. It can be used in the homes of handicapped children in the lounge.

Using the product

Protac MyBaSe® has 2 sides – active and passive. The active side is used when there is a need for vestibular stimulation (balance), and the passive side for when stability is needed.

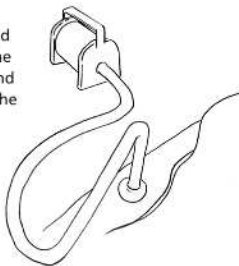
For extra stimulation, a helper can increase movement when the active side of the mattress is used.

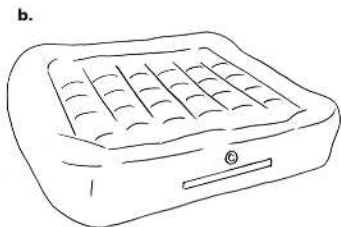
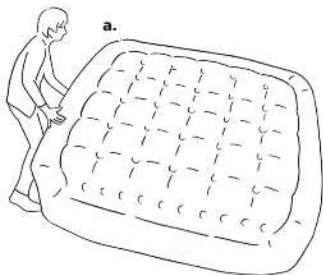
Follow this procedure:

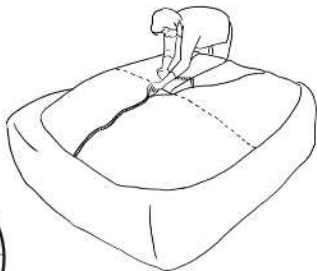
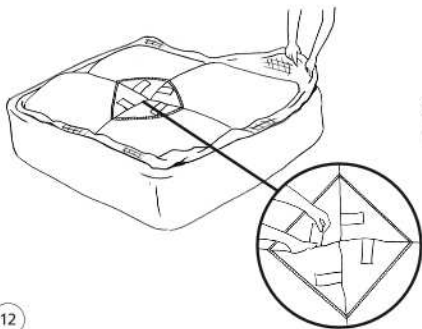
1. When inflating the mattress, open the valve by pressing it inwards and turning it slightly to the right. Attach the threaded end of the hose direct to the mattress, and the other end to the pump. When the mattress is fully inflated, unscrew the hose and seal the valve by pressing it inwards and turning it slightly to the left.

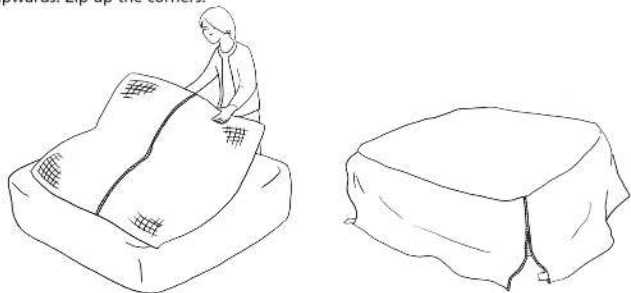
2. Once the mattress is inflated – choose a. the stable side or b. the active side.

3. Place the cover on the mattress with the net side facing upward. Place the 4 ball bags on top with the white zippers facing upwards forming a rectangle in the middle. Join the 4 ball bags using Velcro and place them inside the net cover and zip up.

4. Turn the entire net cover containing the ball bags around. The outer cover is now facing upwards. Zip up the corners.

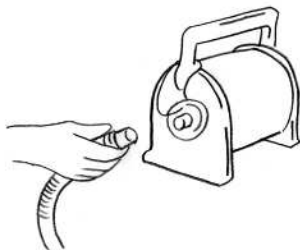
GB

5. There are four handles in the mattress concealed by zippers, making it easy to access the air valve and handles without having to remove the cover.

6. When deflating the mattress, open the valve by pressing it inwards and turning it slightly to the right. Attach the threaded end of the hose to the mattress. Place the nozzle supplied on the pump before attaching the hose.

Accessories

Protac Incontinence Cover

Using the product

The Protac Incontinence Cover is used over the outside of the blanket if the patient is incontinent (involuntary urination).

General precautions - Safety

- ⚠ • The same risk precautions apply as for the use of a plastic bag:
 - A Protac MyBaSe® used with an incontinence cover **must not** be used for children without supervision.
 - A Protac MyBaSe® with incontinence cover **must not** be used for patients unable to care for themselves without supervision.
- ⚠ • **Must not** be used near an open flame.
- ⚠ • **Do not** allow the user to lie with their head down on the mattress.
- ⚠ • You **must not** crawl under the Protac incontinence cover.

Product service life

The product has an expected service life of 5 years.

This is based on correct use, cleaning, washing, service and maintenance.

Products used beyond the aforementioned service life should be assessed for further use.

Responsibility for using the product beyond the stated service life rests with the user.

Cleaning

The air mattress can be cleaned with a cloth wrung in hot water. The cover can be washed. The regularity and need for washing depends on use. See below for washing and cleaning instructions.

Protac MyBaSe®

- The cover can be machine washed in a min. 6 kg capacity washing machine with washing powder with no bleaching agents - see also the washing instructions on the cover.

To remove spots using a spot remover, always try first on an area which is not visible to test for damage or discolouring. Protac does not recommend the use of spot remover.

Protac cannot be held liable for incorrect washing and maintenance.

Maintenance

We recommend washing the cover regularly.

Transport and storage

Protac MyBaSe® must be installed before it is ready for use.

For transport purposes Protac MyBaSe® is packed in 3 carrying bags and protective plastic to protect Protac MyBaSe® against dirt and moisture.

When not in use, Protac MyBaSe® should be stored safely to protect it against dirt and moisture.

When not in use, Protac MyBaSe® can stand against the wall.

Disposal

The Protac MyBaSe® can be broken down into the following components for disposal: inflatable mattress and cover.

All components can be disposed of as combustible refuse.

Always dispose of refuse according to local rules.

Guarantee

2 year guarantee against manufacturing defects.

The guarantee will be invalid in the event of the following:

- over-use
- misuse
- non-maintenance or incorrect maintenance

GB

Traceability: The month and year of manufacturer are marked on the washing instructions label.

Key to symbols for cover

Example only. Always follow the recommended care instructions on the label found inside the product

Washing bowl: Maximum washing temperature (see care label)

Warning:
See under **Safety**

Triangle:
Do not bleach

User information

Square with a circle:
tumble drying possible
(see care label)

For indoor use

CE approval according to Directive concerning medical devices 93/42/EC

Iron:
Do not iron

Manufacturer

Circle:
Do not dry clean

Product serial number (see care label)

Technical specifications

Dimensions:

Protac MyBaSe®:

length: 200 cm.

width: 200 cm.

The air pump is CE-approved.

Standards: Comply with general requirements for aids acc. to the following harmonized standards:

Clinical evaluation: EN 12182:2012

Risk analysis: DS/EN ISO 14971

Service information

No other balls can be used in a Protac MyBaSe® than those used by the manufacturer.

Guide

- to daily inspection of the Protact MyBaSe before use

No.	Check	Remarks
1	<ul style="list-style-type: none">▪ Air mattress intact when inflated	
2	<ul style="list-style-type: none">▪ No sign of holes, damage to the textile or seams	

GB

Protac A/S

Niels Bohrs Vej 31 D, Stilling

DK-8660 Skanderborg

Tel.: 86 19 41 03

e-mail: protac@protac.dk

Dealer list:

Please refer to: www.protac.dk/kontakt

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